




Advanced Injection Belly

LF09900U

INSTRUCTION CARD

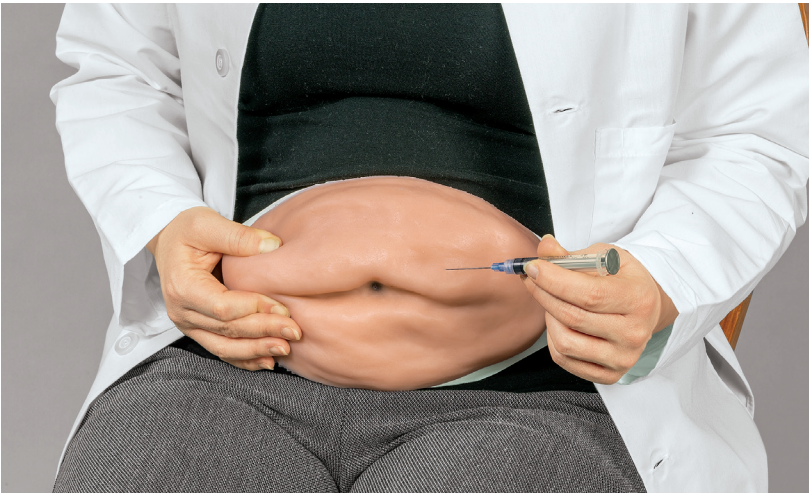

ABOUT THE SIMULATOR

The Advanced Injection Belly is a soft stomach replica developed by Nasco. Made from lifelike skin material, it can be pinched and then injected for a realistic learning experience. Use the straps to apply the belly to a standardized patient or simulator. With little to no maintenance required, this trainer will enhance training scenarios or patient education. The Advanced Injection Belly is made to sustain multiple punctures without showing holes. Protective backing will keep needles from passing through the simulator.

5-Year Warranty

GENERAL INSTRUCTIONS FOR USE

Use for the following types/styles of injection practice:

1. Inserting and rotating the infusion set for insulin pumps
2. Pen-style injections
3. Subcutaneous injections

Note: Nasco recommends air injections only.

CARE AND MAINTENANCE

Clean the Advanced Injection Belly with warm, soapy water. Lightly dry trainer and backing with paper towel. For stubborn dirt, use Nasco Cleaner (LF09919U sold separately). Apply Nasco Cleaner to the soiled area and wipe clean with a clean, dry towel.

A light dusting of baby powder will help maintain the soft, lifelike skin feeling of the simulator and prevent it from becoming tacky.

CAUTIONS

Standard needle handling procedures should be observed when using this trainer. Ink, ballpoint pen, certain printed papers, and plastics will leave an indelible stain when in contact with the simulator.

Actual product may vary slightly from images. Nasco reserves the right to change product color, materials, supplies, or function as needed.

SUPPLIES/REPLACEMENT PARTS

LF09919U Nasco Cleaner

Other Available *Life/form*® Simulators

LF00694U Intramuscular Injection Sites Poster, Chart, and TearPad™ Education Set

LF00952U Injection Training Models, Light

LF00953U Injection Training Models, Dark

LF00961U Intramuscular Injection Simulator

LF01184U Venatech IM & SubQ Simulator, Light

LF01185U Venatech IM & SubQ Simulator, Dark

LF01008U Intradermal Injection Simulator



LF00694U



ONLINE:
E-MAIL:
PHONE:
ADDRESS:

NascoHealthcare.com
lifeform@enasco.com
1.800.431.4310
901 Janesville Ave.
Fort Atkinson, WI 53538